

INTENDED USE

The Allergen (Beef) Rapid Test is a rapid test for the qualitative determination of Beef specific Immunoglobulin E (sIgE) in human serum, plasma or whole blood. The test, in conjunction with other clinical observations, is intended to identify the patient whose allergic symptoms maybe mediated by Beef-specific immunoglobulin E (IgE) Type I hypersensitivity.

INTRODUCTION

Allergy is a common health problem, affecting approximately 20-25% of people with immediate-type hypersensitivity reactions that manifest in the form of rhinitis, urticaria, dermatitis, gastrointestinal illness, wheezing and rarely anaphylactic shock. The term allergy is often used for type I hypersensitivity reactions (immediate type reactions), whose symptoms generally occur within 30-60 minutes after contact with the allergen.

The allergens causing type I hyper sensitivity reactions are mostly proteins derived from the natural environment e.g. plant pollen, animal hair, food, mites, and insect venoms. A characteristic of type I allergies is the involvement of allergen specific immunoglobulins (antibodies) of class E (sIgE). Hence, the detection of sIgE is an important tool of modern allergy diagnostics.

Beef allergy causes an overreaction of the immune system which may lead to severe physical symptoms for millions of people around the world. Beef allergy appears mainly, but not exclusively, in adult.

PRINCIPLE

The Allergen (Beef) Rapid Test Device has been designed to detect Beef sIgE through visual interpretation of color development in the internal strip. The membrane was immobilized with streptavidin on the test region, the conjugate pad was pre-coated with colored anti-IgE antibody colloidal gold conjugates and the sample pad was pre-coated with biotinylated Beef Protein. After specimens were added, the gold-conjugates move along the membrane chromatographically by capillary action and antibodies get to the test region. If sufficient Beef sIgE is present in the sample, it will react with biotinylated Beef protein in sample pad, the mixture then migrates through conjugate pad by capillary action and interact with colored anti-IgE antibody colloidal gold conjugates, form a complex. Then the complex moves to the membrane, and combine with streptavidin. As a result, a colored band will form at the test region of the membrane.

If no Beef sIgE is present in the sample, biotinylated allergen pre-coated on the sample pad will bind to streptavidin immediately, so there is no colored line at the test region of the membrane.

Therefore, the colored band on the test region indicates a positive result. And appearance of a colored band at the control region serves as a procedural control. This indicates that proper volume of specimen has been added and membrane wicking has occurred.

MATERIALS

Materials Provided

- Test Devices
- Package insert
- disposable pipette
- Whole Blood Buffer

Materials Required but Not provided

- Positive and negative controls
- Timer

PRECAUTIONS

- For professional in vitro diagnostic use only.
- Do not use after the expiration date indicated on the package. Do not use the test if the foil pouch or canister is damaged. Do not reuse tests.
- This kit contains products of animal origin. Certified knowledge of the origin and/or sanitary state of the animals does not completely guarantee the absence of transmissible pathogenic agents. It is, therefore, recommended that these products be treated as potentially infectious, and handled by observing usual safety precautions (e.g., do not ingest or inhale).
- Avoid cross-contamination of specimens by using a new specimen collection container for each specimen obtained.
- Read the entire procedure carefully prior to testing.
- Do not eat, drink or smoke in the area where specimens and kits are handled. Handle all specimens as if they contain infectious agents. Observe established precautions against microbiological hazards throughout the procedure and follow standard procedures for the proper disposal of specimens. Wear protective clothing such as laboratory coats, disposable gloves and eye protection when specimens are assayed.
- Humidity and temperature can adversely affect results.
- Used testing materials should be discarded in accordance with local regulations.

STORAGE AND STABILITY

- The kit should be stored at 2-30°C until the expiry date printed on the sealed pouch or canister.
- The test must remain in the sealed pouch or closed canister until use.
- **Do not freeze.**
- Kits should be kept out of direct sunlight.
- Care should be taken to protect the components of the kit from contamination. Do not use if there is

evidence of microbial contamination or precipitation. Biological contamination of dispensing equipment, containers or reagents can lead to false results.

SPECIMEN COLLECTION AND STORAGE

- The Allergen (Beef) Rapid Test Device is intended for use with human whole blood, serum and plasma only.
- Only clear, non-hemolyzed specimens are recommended for use with this test. Serum or plasma should be separated as soon as possible to avoid hemolysis.
- Perform testing immediately after specimen collection. Do not leave specimens at room temperature for prolonged periods. Serum and plasma specimens may be stored at 2-8°C for up to 3 days. For long term storage, specimens should be kept below -20°C. Whole blood collected by venipuncture should be stored at 2-8°C if the test is to be run within 2 days of collection. Do not freeze whole blood specimens. Whole blood collected by finger stick should be tested immediately.
- Bring specimens to room temperature prior to testing. Frozen specimens must be completely thawed and mixed well prior to testing. Avoid repeated freezing and thawing of specimens.
- If specimens are to be shipped, pack them in compliance with all applicable regulations for transportation of etiologial agents.

PROCEDURE

Bring tests, specimens, and/or controls to room temperature (15-30°C) before use.


1. Remove the test from its sealed pouch, and place it on a clean, level surface. Label the test with patient or control identification. For best results, the assay should be performed within one hour.
2. For Serum or Plasma specimens: Hold the dropper vertically and transfer 2 drops of serum or plasma (approximately 50 µL) to the specimen well (S) of the test device, then start the timer..

For Whole Blood specimens: Hold the dropper vertically and transfer 3 drops of whole blood (approximately 75µL) to the specimen well (S) of the test device, then add 1 drop of buffer (approximately 40µL) and start the timer.


Avoid trapping air bubbles in the specimen well (S), and do not add any solution to the result area.

As the test begins to work, color will migrate across the membrane.

3. Read the result visually at 10min. Do not interpret the result after 15 minutes.

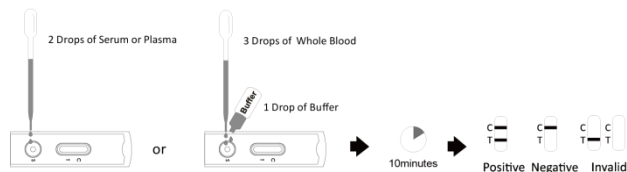

INTERPRETATION OF RESULTS

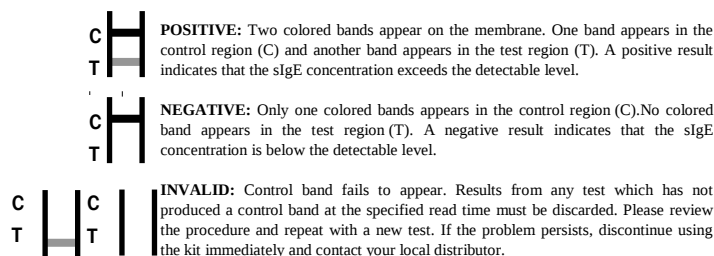

- Note
- 1) The intensity of color in the test region (T) may vary depending on the concentration of analytes present in the specimen. Therefore, any shade of color in the test region (T) should be considered positive. Please note that this is a qualitative test only, and cannot determine the concentration of analytes in the specimen.
 - 2) Insufficient specimen volume, incorrect operating procedure or expired tests are the most likely reasons for control band failure.

QUALITY CONTROL

- Internal procedural controls are included in the test. A colored band appearing in the control region (C) is considered an internal positive procedural control, confirming sufficient specimen volume and correct procedural technique.
- External controls are not supplied with this kit. It is recommended that positive and negative controls be tested as a good laboratory practice to confirm the test procedure and to verify proper

test performance.

LIMITATIONS OF THE TEST

- The Allergen (Beef) Rapid Test Device is for professional in vitro diagnostic use, and should be only used for the qualitative detection of allergen sIgE.
- The Allergen (Beef) Rapid Test Device will only indicate the presence of sIgE in the specimen and should not be used as the sole criteria for the diagnosis of Allergy.
- If the test result is negative and clinical symptoms persist, additional testing using other clinical methods is recommended.
- As with all diagnostic tests, a confirmed diagnosis should only be made by a physician after all clinical and laboratory findings have been evaluated.

PERFORMANCE CHARACTERISTICS

Sensitivity

The analytical sensitivity of The Allergen (Beef) Rapid Test is 0.7IU/mL

Accuracy

A multi-center clinical evaluation was conducted comparing results obtained using The Allergen (Beef) Rapid Test to another commercially available Allergen Rapid Test. The results of the study, which included 62 serum specimens, demonstrated 96.9 % accuracy of The Allergen (Beef) Rapid Test when compared to Allergen Western Blot

Allergen (Beef) Rapid Test vs. Western Blot				
Method	Western Blot			Total results
	Results	Positive	Negative	
	Positive	25	1	26
Allergen (Beef) Rapid Test	Negative	1	38	39
	Total Results	26	39	65

Relative Sensitivity: 96.2% (81.1%-99.3%)

Relative Specificity: 97.4% (86.8%-99.5%)

Overall Agreement: 96.9% (89.5%-99.2%)

Interference Testing

The following substances were added to Beef sIgE free serum and serum spiked with 0.7IU/mL Beef sIgE. None of the substances interfered with the assay at the listed concentrations.

Acetaminophen	20 mg/dL	Acetylsalicylic Acid	20 mg/dL
Ascorbic Acid	20 mg/dL	Atropine	20 mg/dL
Caffeine	20 mg/dL	Gentisic Acid	20 mg/dL
Glucose	2 g/dL	Hemoglobin	1 mg/dL

GLOSSARY OF SYMBOLS

ρ	Catalog number	0	Temperature limitation
ι	Consult instructions for use	Λ	Batch code
⌈	In vitro diagnostic medical device	ε	Use by
μ	Manufacturer	⌈	Contains sufficient for <n> tests
σ	Do not reuse	A	Authorized representative in the European Community
Υ	CE marking according to IVD Medical Devices Directive 98/79/EC		